


Hurricane™ RX

Biliary Balloon Dilatation Catheter

**Boston
Scientific**
Advancing science for life™

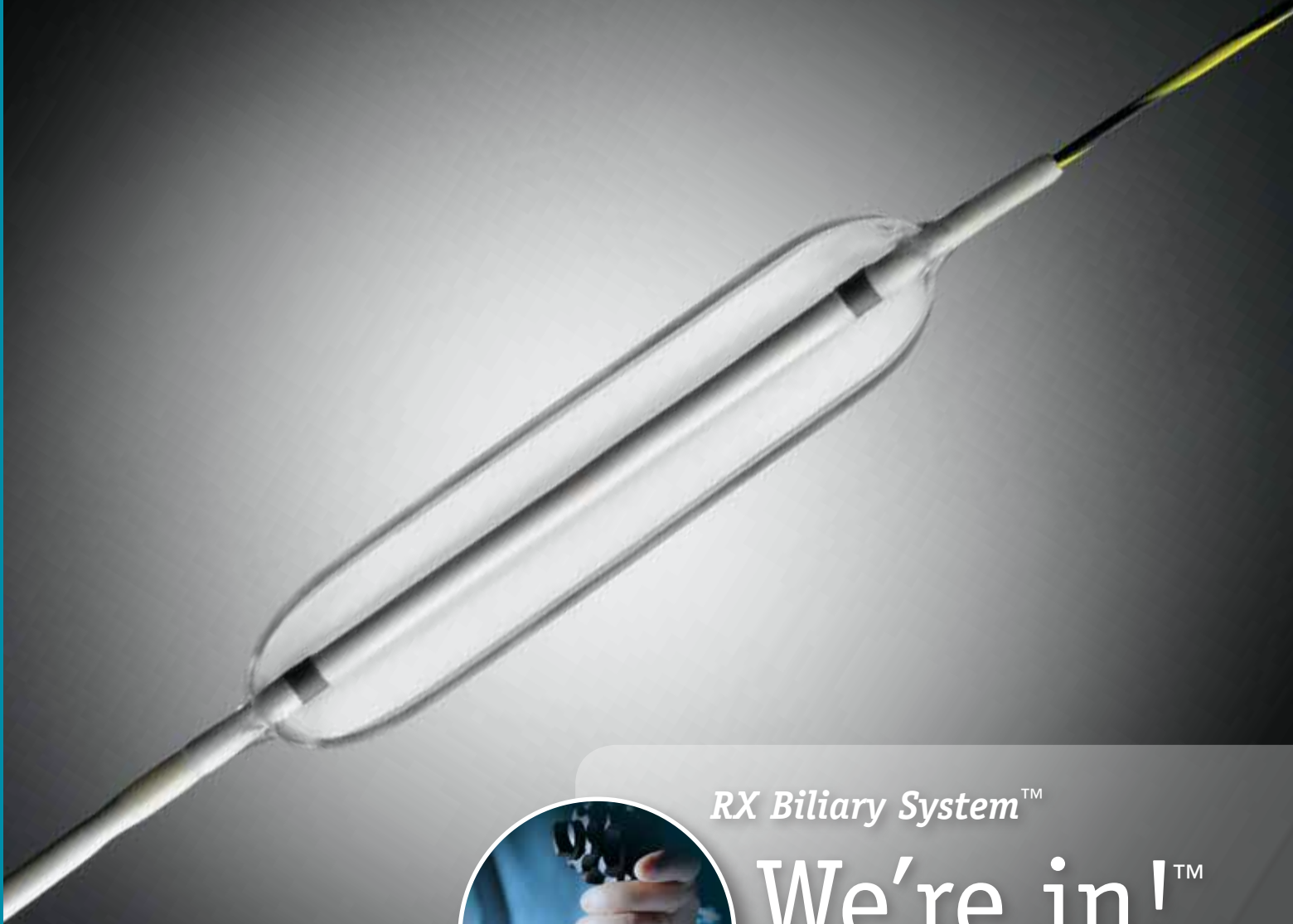

RX Biliary System™

We're in!™

Hurricane™ RX

Biliary Balloon Dilatation Catheter

The Hurricane RX Biliary Balloon Catheter is a double-lumen catheter for endoscopic dilation of strictures of the biliary tree and the Sphincter of Oddi. The balloon is available in a range of sizes to accommodate patient anatomy and stricture size.

RX Biliary System™

We're in!™


The RX Biliary System helps put the physician in control of the guidewire. Using the RX Locking Device, the physician can secure the guidewire to maintain access while exchanging devices.

The RX Biliary System and its compatible devices offer a broad portfolio of products for pancreatico-biliary stone and stricture management.

Excellent Stricture Passage

Low-profile, 5.8Fr balloon shaft is designed to reduce resistance and increase pushability when crossing strictures


- 4Fr tip is designed for less traumatic advancement across the stricture
- Stainless steel stiffening stylet inserted in the proximal portion of the catheter shaft to augment pushability
- RX catheter design allows guidewire locking to secure access

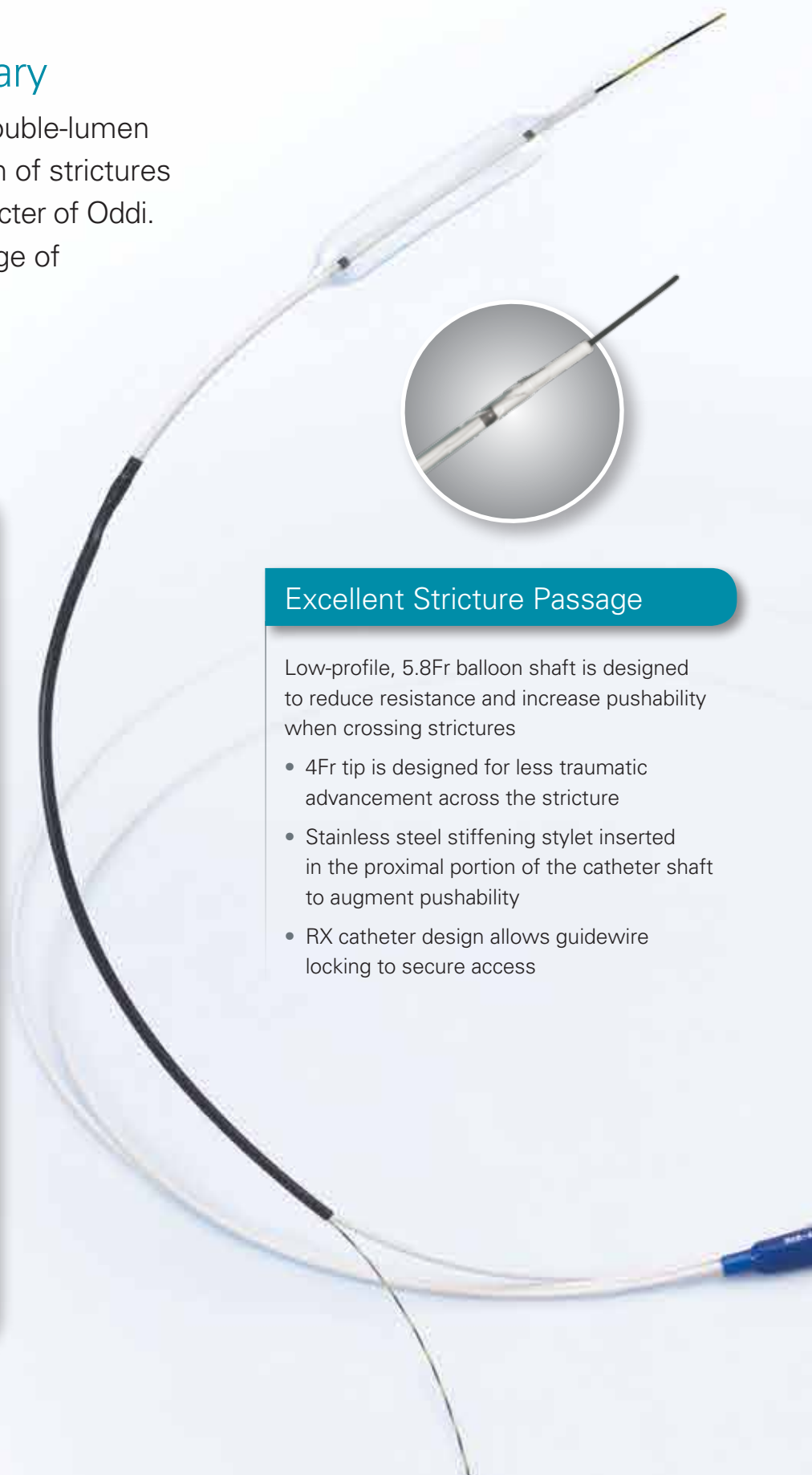

Designed for Accurate Positioning

Hurricane™ RX Balloon offers multiple options for visualization

- Two radiopaque markers aid in balloon positioning
- Contrast may be injected with guidewire in place

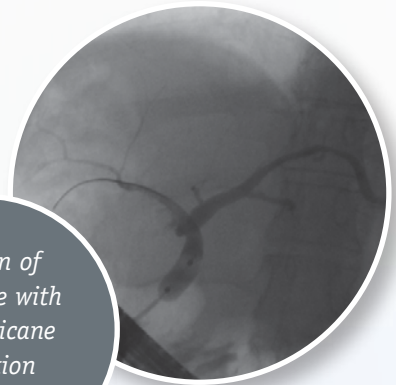

Dilation of a stricture with the Hurricane Dilatation Balloon

Image courtesy of Mario Traina, MD, ISMETT Hospital, Palermo, Italy

Range of Balloon Sizes

Options of balloon diameters and lengths facilitate selection of appropriately sized balloon for patient anatomy and stricture size

- 2cm and 4cm lengths available
- 4mm, 6mm, 8mm and 10mm balloon diameters available

Balloon Achieves Pressure Up to 12atm

- Long 4cm balloon length is designed to offer stability during stricture dilation
- Balloon is recommended for use with the Encore™ 26 Inflation Device which has been developed for rapid deflation to potentially reduce patient discomfort


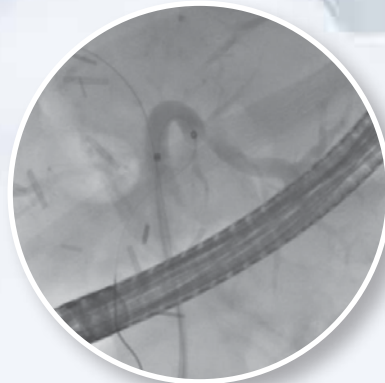


Image courtesy of Sri Komanduri, MD, Feinberg School of Medicine, Northwestern University

Hurricane™ RX

Biliary Balloon Dilatation Catheter

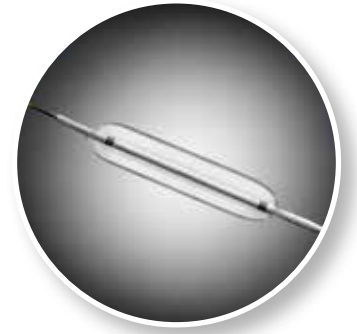

Hurricane RX Biliary Balloon Dilatation Catheter

Order Number	GTIN	Balloon Inflated O.D. (mm)	Balloon Length (cm)	Catheter O.D. (Fr)	Catheter Length (cm)	Maximum Inflation Pressure (atm)	Minimum Channel Size (mm)	Recommended Guidewire (in)	Packaging
M00545890	08714729283805	4	2	5.8	180	12	3.2	.035	Each
M00545900	08714729283812	4	4	5.8	180	12	3.2	.035	Each
M00545910	08714729283829	6	2	5.8	180	11	3.2	.035	Each
M00545920	08714729283836	6	4	5.8	180	11	3.2	.035	Each
M00545930	08714729283843	8	2	5.8	180	11	3.2	.035	Each
M00545940	08714729283850	8	4	5.8	180	11	3.2	.035	Each
M00545950	08714729283867	10	2	5.8	180	8	3.2	.035	Each
M00545960	08714729283874	10	4	5.8	180	8	3.2	.035	Each

*Recommended Guidewire: .035" Hydra Jagwires Guidewire.

Encore™ 26 Inflation Device

Order Number	GTIN	Description	Packaging
M00566670	08714729755241	Encore 26 Inflation Device	Each

Boston Scientific
Advancing science for life™

Boston Scientific Corporation
300 Boston Scientific Way
Marlborough, MA 01752-1234
www.bostonscientific.com

Ordering Information
1.888.272.1001

©2016 Boston Scientific Corporation
or its affiliates. All rights reserved.

ENDO-334301-AB May 2016

All trademarks are the property of their respective owners.

Indications, Contraindications, Warnings and Instructions for Use can be found in the product labeling supplied with each device.

Caution: Federal (U.S.) law restricts this device to sale by or on the order of a physician.